


# emed, pedar, and pliance Recorder software

novel Recorder software option perfectly compliments novel emed, pedar, and pliance pressure measurement software. With the novel Recorder option, digital video is recorded along with pressure data. Both sets of data can be immediately and synchronously played back. This makes Recorder software the perfect tool for display, presentations, education, gait training, and analysis.

The Recorder software and digital video camera setup is effortless. Simply connect a digital video camera(s) to the computer via the FireWire (also known as i.link or IEEE 1394) port, as normally done to record video directly to the PC, set up the novel pressure measurement equipment, and start the novel measurement software. The digital video is automatically displayed on the pressure measurement software screen. With the ability to record pressure along with up to three cameras' video, you are guaranteed not to miss any of the action.

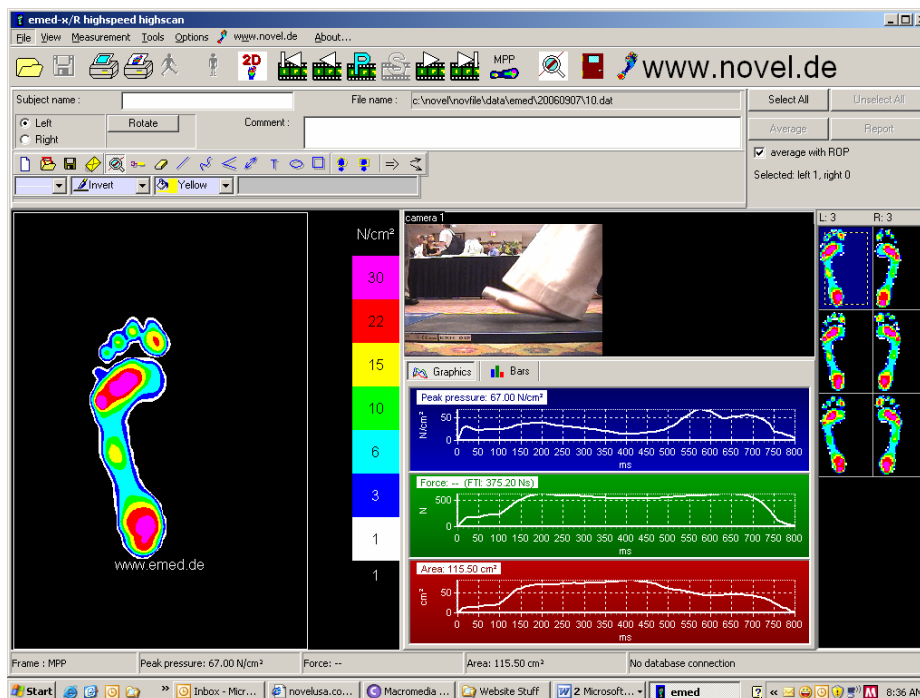

emed/Recorder screen

After data collection, the synchronized recordings can be played back in the measurement software or imported into novel player software to create a media file (.avi) of the pressure and the video. The video from player can be used with PC or Macintosh media players, in presentations, or anywhere standard video AVI files can be viewed.